

#### **Gracey Style Offerings:**

**Standard** - For removal of fine-moderate calculus

- Optimal for pockets 4mm and less
- Greater tactile sensitivity & flexibility than rigid curettes

**Rigid** - For removal of heavy calculus


- Less flexible
- Thicker shank

After Five - Lower terminal shank extended 3mm

- Optimal for pocket depths greater than 4mm
- Blade is 10% thinner to allow for insertion with less tissue distention
- Provides improved access to root surfaces

Mini Five - Ultra slender blades, 50% thinner


- Reduces tissue distention
- Eases subgingival insertion
- Elongated terminal shank for access to pockets & root surfaces greater than 4mm.

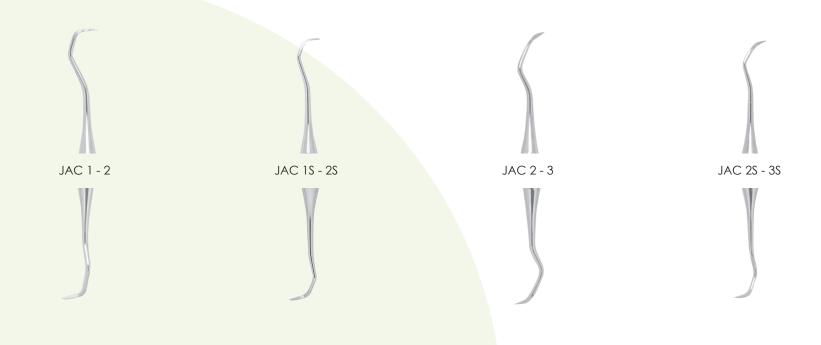

### **PERIODONTAL**

#### **TITANIUM IMPLANT SCALERS/CURETTES**



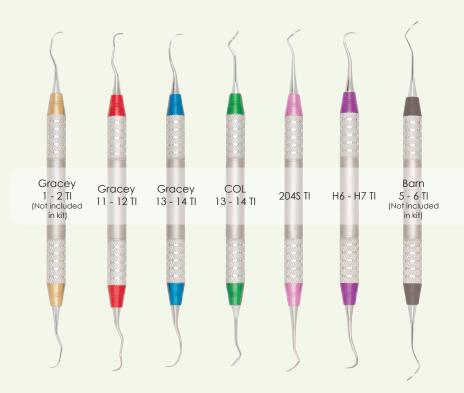

- Color-coded handles
- Made from 6AL-4V titanium alloy
- Fully autoclavable
- Guaranteed sharpness
- Sterilizable titanium tips
- Made in the U.S.A.

#### K1800

Titanium Implant Scaler Kit Includes: Col 13-14 TI 204s TI Gracey 11-12 TI Gracey 13-14 TI H6-H7 TI Sterilization Cassette T-5



## SHARPENATOR

Taking Scaling to New Heights


- Efficiently sharpens instruments without the need of two stones
- Coarse on one side and medium on the other
- The coarse side quickly re-shapes your blade while the medium side delivers a precision edge
- Made from a synthetic material, the stone is autoclavable, non-loading, requiring no oil
- Perfect for chair-side use
- Takes up no counter space

#### The Sharpenator Kit

Includes one of each:
Medium (white)/ Coarse (pink)
Sharpenator Stone, Ergonomic Test
Stick, Sharpenator Chair-side Disposable
Sharpening Guide



Aluminum oxide stones are the most popular sharpening stones in dental schools across the nation. The most preferred stone is the Edgemate<sup>TM</sup> II which offers both coarse and medium sides. All Edgemate<sup>TM</sup> sharpening stones are easy to care for and do not require any lubrication (including water).

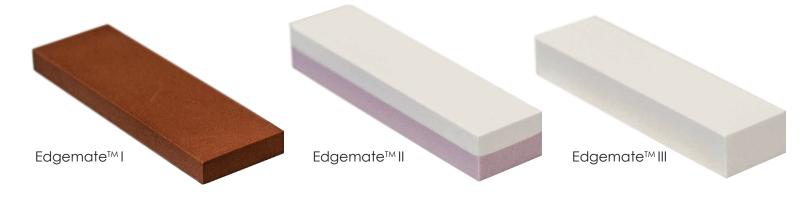


Three grits available:
Fine
Medium
Coarse
Size: 3" x 1" x 1/4"


Coarse on one side and medium on the other. The coarse side quickly re-shapes your blade while the medium side provides a precision edge.


Size: 4" x 1" x 1/2"


Medium grit stone provides a precision edge every time.


Size: 4" x 1" x 1/2"

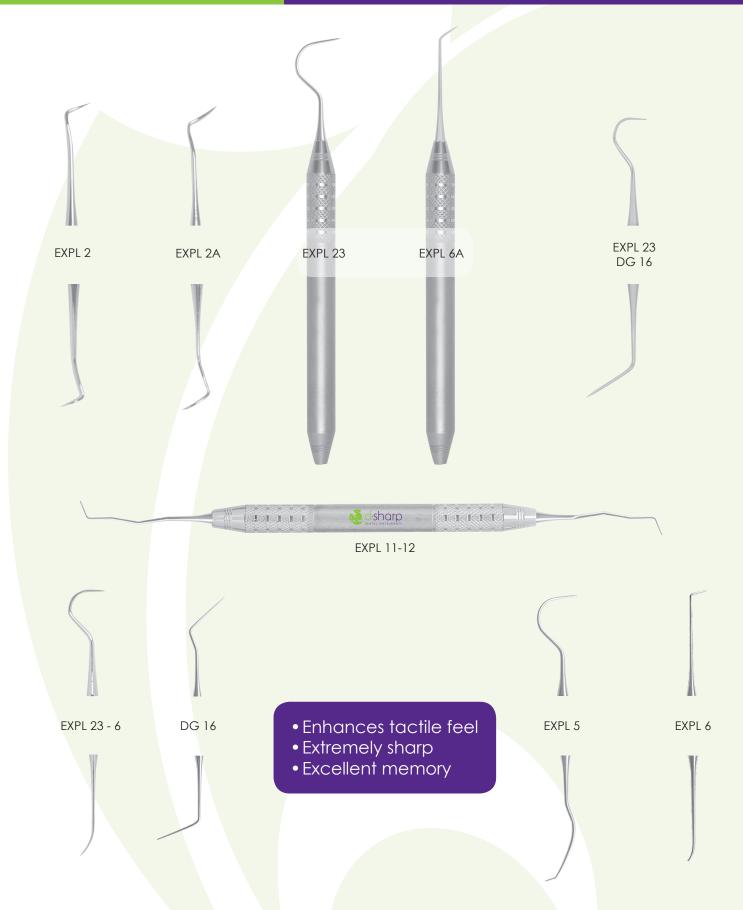

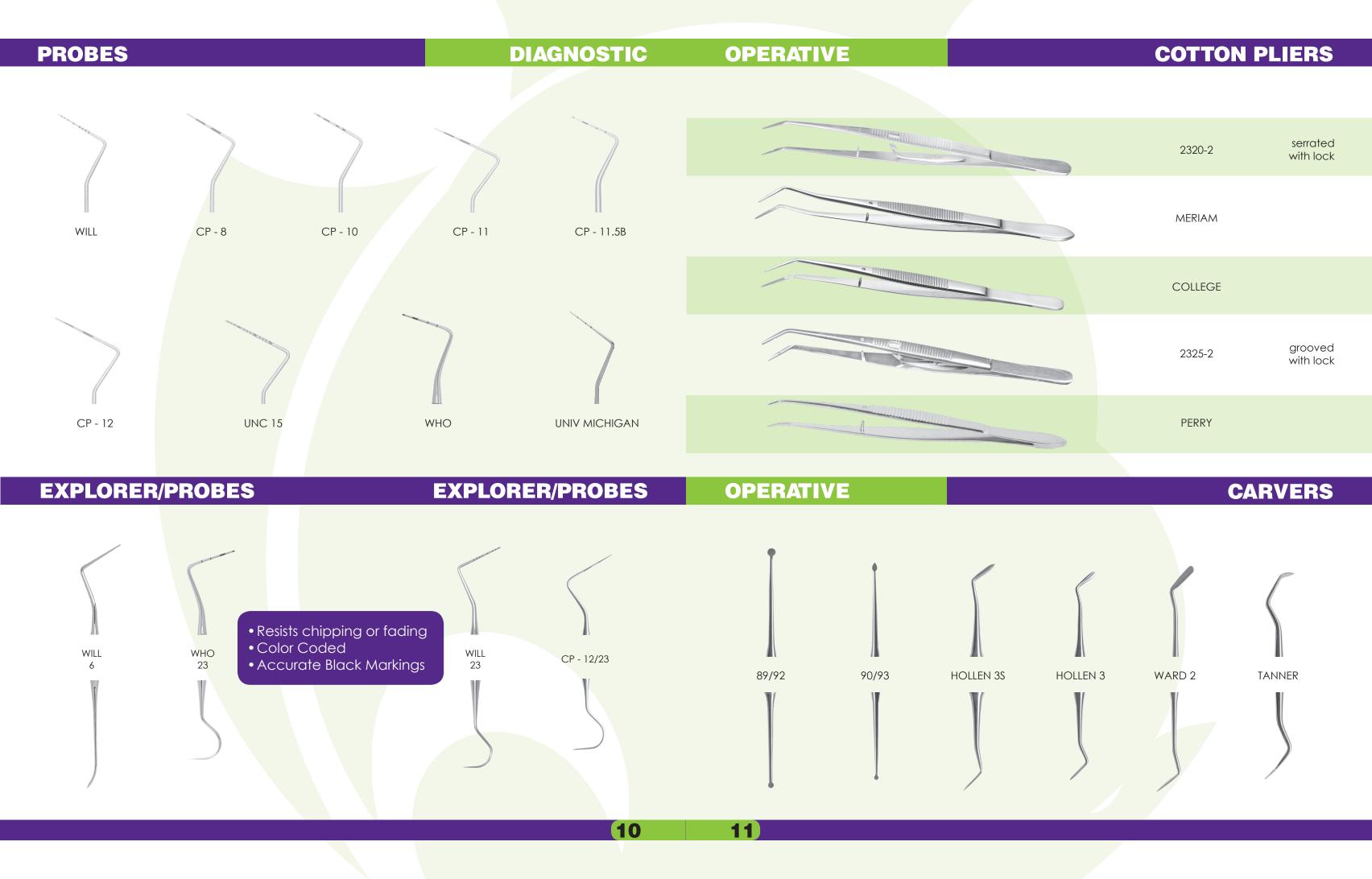

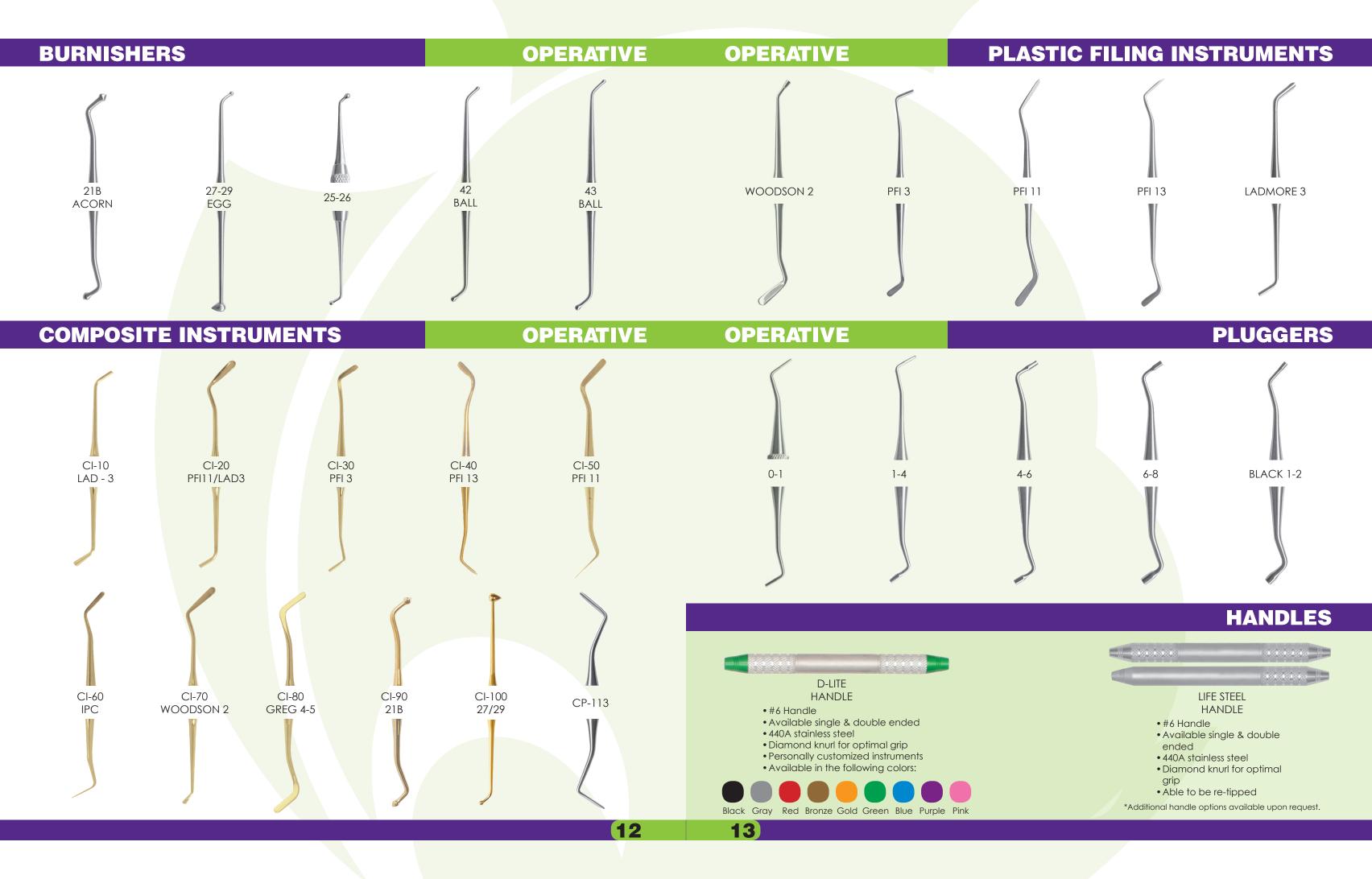

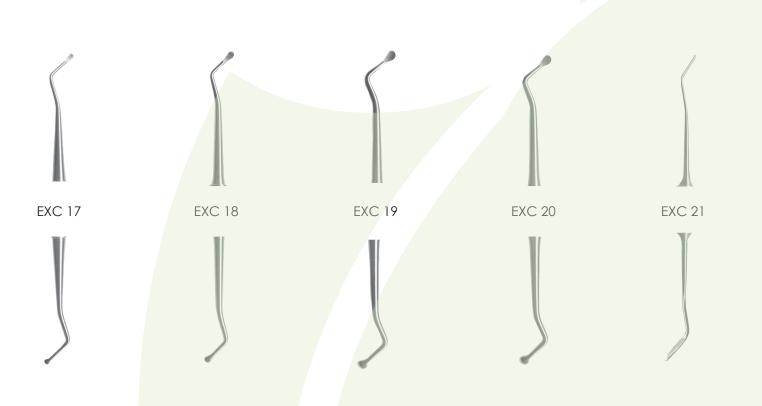

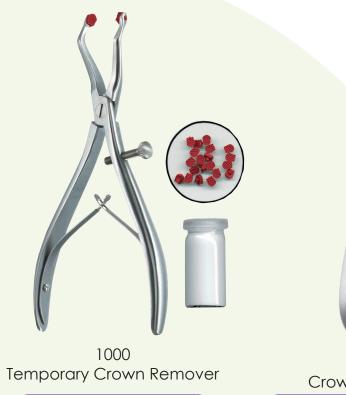

- Spring & screw set prevents crushing of crown
- 20 Rubber Inserts 10g of emery powder




- Spreads the crown & breaks the cement seal
- Unique tip design enables usage throughout entire mouth Round handle provides optimal

maneuverability & leverage



1010 Crown Remover

- Crown spreading forceps
- Spreads the crown & breaks the cement seal
- Built in spring provides even transmission of force
- Eliminates shattering
- Eliminates pressure in patient's mouth

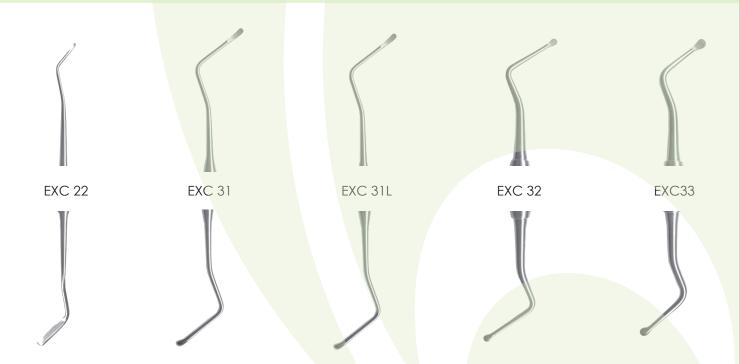

### **STERILIZATION CASSETTES**





T - 5

Custom inserts are made for all sterilization cassettes to hold instruments securely and maximize sterilization.

14

# Blancat Just one hour away from a brighter, healthier smile



## Congratulations!

You are one step closer to getting the healthier and happier looking smile you deserve.

Introducing the BLANCA Tooth Whitening System designed for people with an active life style who want quick, successful whitening results. Our convenient tooth whitening components were developed with you in mind, for simply whiter teeth! Try one of our products today and see why so many people are enjoying brighter, healthier smiles in a fraction of the time.

#### **Professional Strength Whitening**

Our whitening light utilizes the most advanced LED light technology to give you brilliant, professional results.

- Achieve maximum whitening results with two 20 minute sessions
- Remove stains that have built up over time
- State of the art LED technology activates gel with no heat, no UV light, and no tooth dehydration
- Whiten your teeth 3 5 shades

#### **Naturally White**

Our gel is a proprietary blend of Hydrogen Peroxide and Carbamide Peroxide with a neutral PH that will accelerate whitening and not harm your enamel.

- Pre-filled whitening tray
- Foam absorption strips comfortably fit against the upper and lower teeth
- A fraction of the time and cost of other in-office whitening systems and take home trays
- Vacuum-sealed packaging ensures freshness
- Features a special vent for improved breathing and comfort



This take-home Whitening Paste reinforces whitening results each time you apply it to your teeth.

- A naturally derived paste to reduce sensitivity after whitening procedures
- Restores minerals lost during the bleaching process and every day wear
- Contains 8 grams of paste for 30 applications

# BRIGHTEN YOUR TEETH the go

This is the most convenient solution in teeth whitening. Brighten your teeth on-the-go without messy strips, trays or rinsing.

- Easy brush-on applicator pen
- Our fast action whitening formula contains 15% hydrogen peroxide. The most powerful and effective paint on whitening gel available.
- Instantly dissolves film from red wine and coffee
- The 4.5ml, sleek pocket-size "pen" gives you over 30 applications
- Pleasant Mint flavor enhances whitening experience and limits salivation





1. Brush & floss prior to application



2. Place lip guard in mouth



3. Twist pen to activate gel onto brush



4. Paint a thin layer to front surface of teeth

#### **INDEX**

Blanca Teeth Whitening Professional Whitening	16	Plastic F
Tooth Whitening Paste		PFI11/LAD3
Tooth Whitening Pen		PFI3
Burnishers		PFI 13
21B ACORN	10	PFI 11
27-29 EGG		WOODSON 2
25-26		GREG 4-5 21 B
12 BALL		27/29
13 BALL		CP-113
Carvers		PFI 3
39/92	11	PFI 11
90/93		PFI 13
HOLLEN 3S		LADMORE 3
HOLLEN 3	11	<b>Pluggers</b>
VARD 2		0-1
ANNER	11	1-4
Composite Instruments		4-6
CI-10 LAD - 3	12	6-8
CI-20 PFI11/LAD3		BLACK 1-2
CI-30 PFI 3		Probes
CI-40 PFI 13		WILL
CI-50 PFI 11		CP-8
CI-60 IPC		CP-10
CI-70 WOODSON 2 CI-80 GREG 4-5		CP-11
CI-90 21B		CP-11.5B
CI-100 27/29		CP-12
CP-113		UNC 15 WHO
Cotton Pliers		UNIV MICHIGAN
320-2 Serrated with lock	11	Scalers 8
MERIAM		
COLLEGE		204 204s
2325-2 Grooved with lock		204SD
PERRY	11	N 135
Crown & Bridge		N 137
000 Temporary Crown Remover	11	U 15-30
008 Crown Distender	11	H5
010 Crown Remover		H5-33
Excavators		H6-H7
XC 17	14	KB 1
XC 18		MTJ NEVI 1
XC 19		NEVI 1
XC 20	14	NEVI 3
XC 21		NEVI 4
XC 22		GF-5
XC 31		GF-6
EXC 31L		LANGER 1-2
EXC 32		LANGER 3-4
EXC 33	14	LANGER 5-6
Explorers/Probes		LANGER 17-18.
XPL 2		U15
XPL 2A		MC 13S-14S MC 13-14
EXPL 23		MC 17S-18s
XPL 6A XPL 23 DG 16		MC 17-18
EXPL 23 DG 16		BH 1-2
DG 16		BH 5-6
XPL 5		COL 13-14
XPL 6		COL 4R-4L
VILL 6		COL 2R-2L
VHO 23	10	YG 7-8
VILL 23		GR 1-2
CP-12/23	10	GR 3-4
-landles		GR 5-6 GR 7-8
D-LITE HANDLE		GR 9-10
IFE STEEL HANDLE		GR 11-12

Plastic Filling Instruments	Scal
LAD-312	GR 11
PFI11/LAD312	GR 12
PFI312	GR 13
PFI 1312	GR 15 GR 17
PFI 1112	SUB 0
WOODSON 212	SUB 1
GREG 4-512	SUB 1
21 B12	JAC 1
27/2912	JAC 1
CP-11312	JAC 2
PFI 3	JAC 2
PFI 11	Sha
PFI 13	
ADMORE 3	Edgen
Pluggers	Edgen
0-113	Edgen
1-4	Ste
4-6	T-20
5-8	T-16
BLACK 1-213	T-10
Probes	T-5
WILL10	Tita
CP-810	GRAC
CP-1010	GRAC
CP-1110	GRAC
CP-11.5B10	COL 1
CP-1210	204S
JNC 1510	BARN
WHO10	H6-H7
UNIV MICHIGAN10	
Scalers & Curettes	

Scalers & Curettes Cont.	
GR 11-146	
GR 12-13	
GR 13-146	
GR 15-166	
GR 17-18	
SUB 0 6	
SUB 11-126	
SUB 13-14	
JAC 1-2	
JAC 1S-2S	
JAC 2-37	
JAC 2S-3S 7	
Sharpenator	
EdgemateTM I8	
EdgemateTM II	
EdgemateTM III8	
Sterilization Cassettes	
T-20	
T-16	
T-1015	
T-515	
Titanium Implant Scaler/Implants	
GRACEY 1-2 Tl7	
GRACEY 11-12 TI	
GRACEY 13-14 TI	
COL 13-14 TI	
2014 TI 7	



18

#### **INSTRUMENT MAINTENANCE SOLUTIONS**

D-Sharp offers unique instrument life extension programs, new instrument acquisition and management strategies, hands on continuing education for dental hygienists, and other operatory cost containment services to help you reduce operatory costs and increase practice revenue.

**OIM Program:** An instrument management program ensuring that each dental hygienist within your practice will always be using properly shaped and sharp instrumentation. Every 4 months receive 5 sets of instrumentation. Simply return the used sets, using the provided UPS label. Maintenance will be executed on the returned sets and in 4 months time the next shipment will arrive. For \$129/month (\$99/month for multiple operatories) OIM takes the hassle out of instrument management. Call to request information regarding program details and how to sign up.

**VitaSharp** is a sharpening service that will help extend the life of your instruments. With the use of the instruments, the edges become beveled and require re-contouring. Utilizing professional sharpening will extend the life of your instrument and allow you to work efficiently and effectively with instruments that perform like new for an extended period of time and help to provide optimal patient care. (VitaSharp pricing: \$3/per end, allow 3 weeks for service).

**IntroCusp/Re-tipping Service:** State of the art manufacturing process ensures brand new tips are inserted into the existing handles you have already purchased. This option offers a significant cost savings, up to 70%, when compared with a new purchase of scalers or curettes. (\$5.99/end, allow 3 weeks for service)

**Trade In:** By trading in old instruments for new, your practice will have a significant cost savings while decreasing waste and disposal from your practice. When instruments can be refurbished, D-Sharp will renew and donate to those in need. Send in old instrumentation you don't wish to use (min 12) and receive \$6 off each scaler or curette.

#### **ORDERING**

To place an order or inquire about D-Sharp instruments or services please call 877.837.4277. You may also email us at info@dsharp.ca or visit our website www.dsharp.ca

#### **WARRANTY**

All D-Sharp products are warranted to be free from defects in workmanship and materials. Any product which proves to be defective will be repaired or replaced free of charge.\* Warranty doesn't apply to normal wear and tear and for use of instruments beyond limitations.

\*Warranty doesn't cover shipping and handling



877.837.4277

www.dsharp.ca

